


for MH8001 Series

Operation for return a factory setting:


1. Remove battery and AC/DC Power Supply and make sure remove **ALL** the electrodes from skin.
2. Waiting one minute.
3. To insert the battery into the compartment and make sure that you are installing the battery properly.

(make sure remove **ALL** the electrodes from skin)

4. **First** press and hold “CH4 ▼” intensity decrease and “**P-**” buttons at the same time,

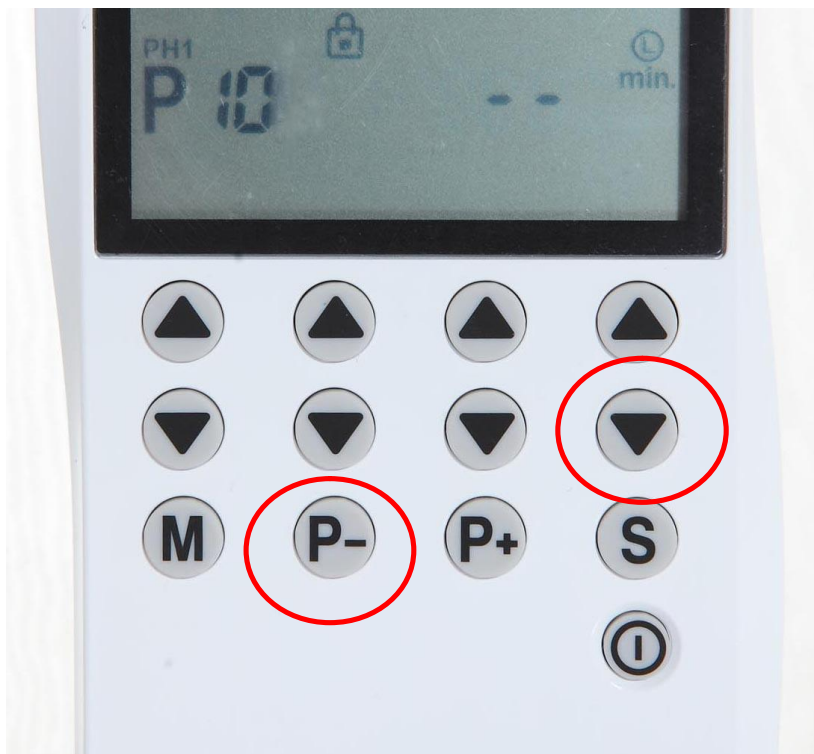

**Then** press and hold power ON/OFF button for 1 second for power on.

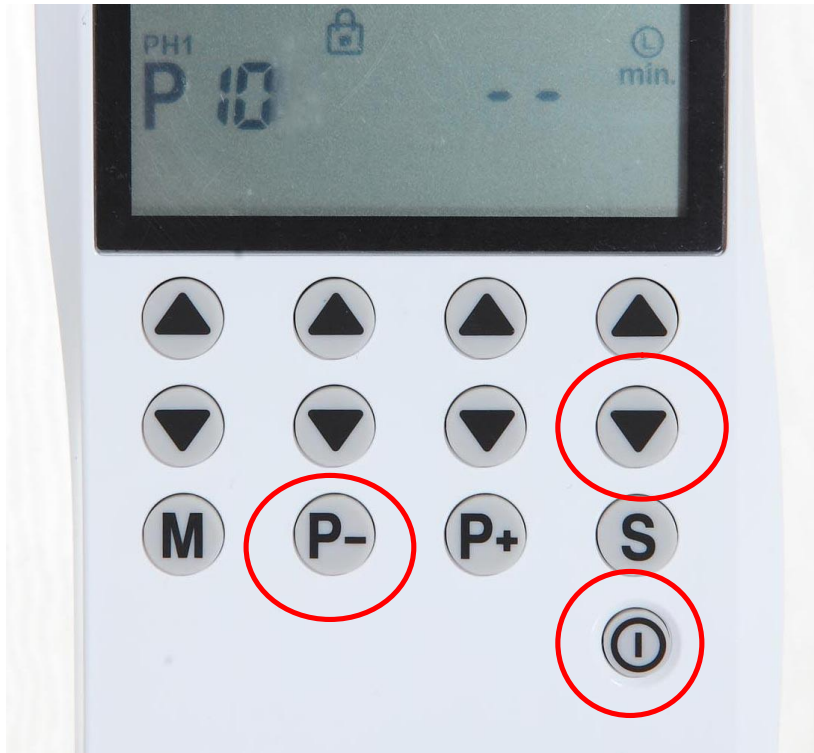

5. The unit will be signaling "Bi" sound about 2 seconds.
6. Please release all buttons.
7. Press the Power ON/OFF button to switch off unit.
8. The unit already returned a factory setting